



μBeam Pro

NEW

The Analyzer for Microscopic beams

PRELIMINARY



The μBeam Pro is a high-precision beam diagnostic system designed for real-time measurement and visualization of micronic continuous wave (CW) or pulsed lasers. It effectively measures beam profiles from fiber optics, laser diodes, and focused infrared (IR) laser beams.

Key Features:

- Capable of measuring beam sizes smaller than 0.5 μm (FWHM).
- Supports both CW and pulsed lasers with an extended working distance.
- Equipped with a specialized IR camera, including sensitivity to 1550 nm, featuring a global shutter and a 1600x1100 pixel resolution.
- The system includes a 100x magnification objective, covering a wavelength range from 480 nm to 1550 nm.
- A unique filter arrangement allows the system to handle higher laser power levels, adding flexibility to the objective lens.



Specifications

Design Wavelength DWL (nm):	480 - 1550
Objective – Infinity Corrected	X100 NIR infinity with tube lens and optical assembly
Numerical Aperture NA:	0.50
Working Distance	12 mm
Attenuation	Built-in thread filter adaptors at the objective's back sight
Interface	USB 3.0 – Windows Operating system
Software Package	Installable detailed GUI for system operation

Camera Sensor Size	14.5 x 9.9 mm
Sensor Resolution	1600 x 1100, Global Shutter
Gain Control	User Selectable
Exposure Speed	32 μsec to 2 sec, user selectable
ADC Dynamic Range	16 bit
Software Features	
Beam Profiling	Complete XY Beam Analysis including 3-D display and sub-micron results
Data Saving	Full data saving and retrieving including video mode
Image Analysis	Sophisticated image analysis and auto contrast stretch
Multiple Beam Analysis	Available

Working Distance

12

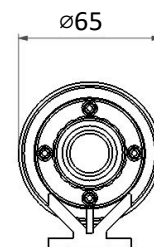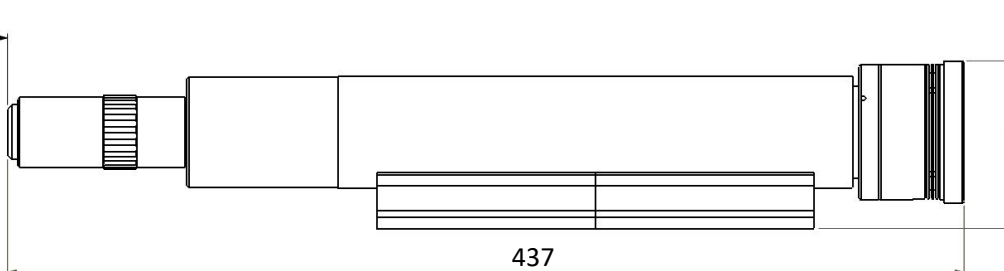

DUMA OPTRONICS LTD.

Website: <http://www.dumaoptronics.com>

E-mail: sales@duma.co.il

October 2024

